GB

Models covered by this information

Product Name	Part No	Standards
Knee Brace C-Frame Pro Carbon #S/M Pair	5017010100	EN 1621-1:2012
Knee Brace C-Frame Pro Carbon #L/XL Pair	5017010101	
Knee Brace C-Frame Pro Carbon #XXL Pair	5017010102	
Knee Brace C-Frame Pro Carbon #S/M Left	5017010110	
Knee Brace C-Frame Pro Carbon #L/XL Left	5017010111	
Knee Brace C-Frame Pro Carbon #XXL Left	5017010112	
Knee Brace C-Frame Pro Carbon #S/M Right	5017010120	
Knee Brace C-Frame Pro Carbon #L/XL Right	5017010121	
Knee Brace C-Frame Pro Carbon #XXL Right	5017010122	

Testing and Approval

Leatt® Knee Brace C-Frame Carbon is a Knee Brace Class I Medical Device incorporating Personal Protective Equipment (PPE) in the form of a knee protector intended for use by motorcyclists' and downhill cyclist. Leatt® Knee Brace C-Frame Carbon's knee protector have been tested against the manufacturer's specifications and EC Type-Examination certification has been issued by Eurofins-Modulo Uno SpA, Torino, Italy; Notified Body No: 0477. The protectors are CE marked to denote compliance with the EU Personal Protective Equipment Directive 89/686/EEC. The Knee Brace has been assessed against the requirements of the following technical specifications:

- EN 1621-1:2012 – "Motorcyclists' protective clothing against mechanical impact - Part 1: Motorcyclists' limb joint impact protectors."
- EN 13688: 2013 – "Protective clothing – General Requirements"
- Medical Device Directive 92-42-EEC of 1993.
- PPE Directive 89/686/EEC

Protection provided

Protector CE certified to PPE Directive 89/686/EEC and Medical Device Directive 92-42-EEC. Knee Protection conforms to Type A Level 1 requirements of EN 1621-1:2012

Leatt® Knee Brace C-Frame Carbon has been tested to EN1621-1:2012 and recorded a mean maximum transmitted force below 35kN from a 50 Joule impact (Level 1).


Use

The Leatt® Knee Brace C-Frame Carbon provides stability to the knee and is ideal for riders taking part in off-road motorcycling, motocross and downhill mtb. It is also CE certified as a knee guard over a defined range of riding positions. Throughout the range the focus is on knee support with a carbon based structure and 3DF foam for knee impact protection. Products meeting the requirements of the technical specifications listed above are intended to reduce the severity or incidence of injuries related to Hyperextension such as Posterior Cruciate Ligament ruptures (PCL) Anterior Cruciate Ligament ruptures (ACL) as well as Valgus related injuries such as medial collateral ligament ruptures and patellofemoral pain. In addition to the aforementioned, the product is intended to reduce the severity or incidence of soft tissue injuries such as bruising, contusions, cuts and abrasions during off-road motorcycling and biking activities, such as Moto-X, Enduro, Downhill MTB, BMX and Adventure Sports.

- Products meeting the requirements of the technical specifications may reduce the severity and/or incidence of some more serious injuries such as minor fractures, ligament ruptures and impact related knee pain.
- Products meeting the requirements of the technical specifications cannot significantly mitigate against major injuries, such as serious fractures resulting from extreme impact forces or torsion/flexion/bending/crushing injuries. Currently available materials and technology do not allow this level of protection to be attained in products which would be acceptable to wear and which would allow physical activities to be conducted and completed satisfactorily

Protection dimensions

The following illustrations demonstrate the dimensions of the zones of protection relative to the coverage provided by the knee protector.

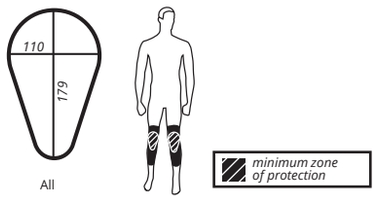Size range

Leatt® Knee braces are available in the following sizes [cm]:
Use the following pictogram to identify the measurements provided.

Upper Leg circumference

Size	cm	inches
S/M	45 - 50	17.5" - 19.5"
L/XL	50 - 58	19.5" - 23"
XXL	58 - 65	23" - 25"

* Measure around the mid part of upper leg

Positioning of Leatt® protection range

No single size of knee brace can accommodate all body dimensions because of variation in height and shape. In addition, the Leatt® Knee Brace C-Frame Carbon can only provide the maximum achievable level of protection when it is correctly positioned and adequately secured on the wearer. Where fitted the C-Frame should always be of a snug fit to the leg. In particular users should avoid choosing too big a protector because it can result in a dangerous riding condition. If the C-Frame does not correctly fit and cover the following areas, an alternative size may be required:

- C-Frame hinge positioned in line with the centre of the knee on the side (in line with the "hinge" of the knee).
- Straps should be secured in the following order to ensure optimal fit:
 1. Calf strap
 2. Lower thigh strap
 3. Upper thigh strap
 4. Bottom calf strap

Adjustment

Unfasten the adjustment straps. Put on the brace and place in position. Pull the straps and adjust to a snug fit. Once the tension of the straps is adjusted, clips can be used to attach and release the straps during subsequent uses. When correctly-adjusted and worn the C-Frame should allow the full range of normal leg movements. With your knee C-Frame in place, check that you are able to safely operate your vehicle's controls. Please ask the retailer or contact Leatt® Corporation if you require further advice on sizing, adjustment or fitting.

Marking and labeling

The following diagrams provide important information about this product and the protection it provides.


CE	The "CE " mark denotes compliance with the health and safety requirements of Annex II of the European Personal Protective Equipment Directive 89/686/EEC.
	The recommendation to read the manufacturer's instructions contained in this booklet
K TYPE A 1 EN 1621-1:2012	Denotes that the product is or incorporates a knee protector ("K - Type A") which meets level 1 of European Standard EN 1621-1:2012 The standard provides for two sizes of protector Type A and Type B representing minimum zones of protection and two levels of protection (impact attenuation), Level 1 and Level 2. Type B protector dimensions are intended to cover the needs of most motorcycle riders. However, for reasons of size, fit, ergonomics or comfort alternative (smaller) Type A protector dimensions may be used. Level of protection is determined by measuring the impact force transmitted in three defined strike areas of the minimum zone of protection, A B and C, after being subjected to a 50 Joule impact under ambient and wet conditions (optionally also hot and/or cold). Level 1 requires a mean impact force transmission of ≤ 35kN, all single strike forces to be ≤ 35kN in strike area A and ≤ 50kN in strike areas B and C. Level 2 requires a mean impact force transmission of ≤ 20kN,all single strike forces to be ≤ 20kN in strike area A and ≤ 30kN in strike areas B and C.

Cleaning

Always follow the washing instructions printed on the label attached to the garments, since when alternative fabrics with specific performance characteristics are used, washing instructions may differ. Touch-and-close fastenings, where fitted, can be cleared of debris using a nail brush. If the knee guard components have become soiled in use, wipe the surface with a damp cloth. Use a soft brush, damped in water, to clean the perforations in the components. Only use tap water. Never use cleaners or detergents, or launder the knee brace, the product will be damaged by the cleaning process and chemicals and the level of protection provided severely reduced.

Common substances applied to this product (solvent, cleaners, hair tonics, etc.) can cause damage that may be invisible to the user and compromise the effectiveness/safety of your Leatt® product. If in any doubt, please contact the manufacturer before cleaning or laundering your protector. Always refer to the material content label on the Leatt® product.

None of the materials used in the construction of the product are known to contain anything that is harmful or might cause an allergic reaction or be dangerous to health.

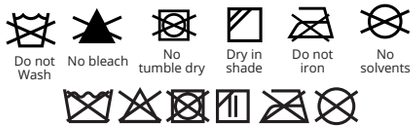clevercare.info

Storage and transportation

Place the Leatt® product on a clean, dry surface and store in a dry, well-ventilated environment away from direct sunlight and extremes of temperature. Do not place heavy objects on top of Leatt® product as this can damage the fabrics and plastic materials. The protection provided by your Leatt® product should not be affected by non-extreme ambient temperature and weather conditions

Maintenance and obsolescence

For maximum performance, the Leatt® product must be inspected prior to each use. Stop using the protector if any components show signs of damage or are worn, cracked, and/or deformed, or if the inside padding has deteriorated. Leatt® Corporation recommends that the Leatt® product is replaced every three (3) years under normal use conditions. Should you experience a collision or if the Leatt® product or any component shows any signs of damage, cracking, deformity and/or the inside padding is deteriorating, return the product to the seller for inspection or destroy and replace it. Please dispose of it in an environmentally-friendly manner. Do not drill holes, paint, cut, heat, bend, apply decals or stickers, or make any other modifications to the knee brace. Doing so may compromise the effectiveness and safety of your Leatt® product.

Use with other Personal Protective Equipment

Leatt® protectors form part of a protective equipment ensemble available for use by participants in off road motorcycling and downhill biking activities. Additional PPE, which it is recommended should be worn during off road motorcycling and downhill mtb activities, includes clothing conforming to EN 13595-1:2002, where applicable to the riding environment (may not be suitable for off-road motorcycling or downhill mtb) or other suitable clothing providing full coverage to the body, arms and legs; impact protectors conforming to EN 1621-1:2012 (where not already a component of your Leatt® protector or other motorcycle clothing); gloves conforming to EN 13594:2002; boots conforming to EN 13634:2010 Always perform any motor sport activity with caution and use/wear appropriate safety equipment.

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Copyright © 2025, Leatt® Corporation
Leatt® Corporation HQ
Leatt® Corporation USA