GB

Models covered by this information

Product Name	Part No	Standards
Knee Brace X-Frame S Pair	5018010101	EN 1621-1:2012
Knee Brace X-Frame M Pair	5018010102	
Knee Brace X-Frame L Pair	5018010103	
Knee Brace X-Frame XL Pair	5018010104	
Knee Brace X-Frame XXL Pair	5018010105	
Knee Brace X-Frame S Left	5018010121	
Knee Brace X-Frame M Left	5018010122	
Knee Brace X-Frame L Left	5018010123	
Knee Brace X-Frame XL Left	5018010124	
Knee Brace X-Frame XXL Left	5018010125	
Knee Brace X-Frame S Right	5018010141	
Knee Brace X-Frame M Right	5018010142	
Knee Brace X-Frame L Right	5018010143	
Knee Brace X-Frame XL Right	5018010144	
Knee Brace X-Frame XXL Right	5018010145	

Testing and Approval

Leatt® Knee Brace X-Frame is a Class I medical device self-certified by Leatt Corporation in terms of Directive 92/42/EEC and registered as EU CA014741 and USA FDA 10048761. The integrated PPE knee protector has been tested against the manufacturer's specifications at test lab Omega S.r.l of Strada Comunale Savonesa 9 – 15050 Tortona – Italy and EC Type-Examination certification has been issued by ITS Testing Services (UK) Limited, Centre Court, Meridian Business Park, Leicester, LE19 1WD, United Kingdom, Notified body No. 0362. The knee protector components are CE marked to denote compliance with the Personal Protective Equipment Regulation (EU) 2016/425. The following European Standards have been used in the evaluation of the knee protector:

• EN 1621-1:2012 – "Motorcyclists' protective clothing against mechanical impact - Part 1: Motorcyclists' limb joint impact protectors."
• EN 13688:2013 – "Protective clothing - General requirements"

Note: To view the Declaration of Conformity to PPE Regulation (EU) 2016/425 please navigate to leatt.com and select the DOCUMENTS tab on your product's page.

Protection provided

Protector CE certified to PPE Directive 89/686/EEC and Medical Device Directive 92-42-EEC. Knee Protection conforms to Type A Level 1 requirements of EN 1621-1:2012

Leatt® Knee Brace X-Frame has been tested to EN1621-1:2012 and recorded a mean maximum transmitted force below 35kN from a 50 Joule impact (Level 1).

Use


The Leatt® Knee Brace X-Frame provides stability to the knee and is ideal for riders taking part in off-road motorcycling, motocross and downhill mtb. It is also CE certified as a knee guard over a defined range of riding positions. Throughout the range the focus is on knee support with a based structure and 3DF foam for knee impact protection.

Products meeting the requirements of the technical specifications listed above are intended to reduce the severity or incidence of injuries related to Hyperextension such as Posterior Cruciate Ligament ruptures (PCL) Anterior Cruciate Ligament ruptures (ACL) as well as Valgus related injuries such as medial collateral ligament ruptures and patellofemoral pain. In addition to the aforementioned, the product is intended to reduce the severity or incidence of soft tissue injuries such as bruising, contusions, cuts and abrasions during off-road motorcycling and biking activities, such as Moto-X, Enduro, Downhill MTB, BMX and Adventure Sports.

- Products meeting the requirements of the technical specifications may reduce the severity and/or incidence of some more serious injuries such as minor fractures, ligament ruptures and impact related knee pain.
- Products meeting the requirements of the technical specifications cannot significantly mitigate against major injuries, such as serious fractures resulting from extreme impact forces or torsion/flexion/bending/crushing injuries. Currently available materials and technology do not allow this level of protection to be attained in products which would be acceptable to wear and which would allow physical activities to be conducted and completed satisfactorily

Protection dimensions

The following illustrations demonstrate the dimensions of the zones of protection relative to the coverage provided by the knee protector.

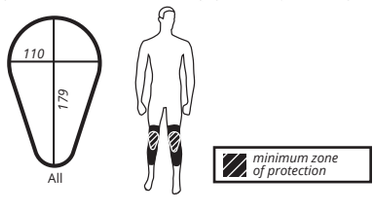Size range

Leatt® Knee braces are available in the following sizes [cm]:

Use the following pictogram to identify the measurements provided.

Knee width

Size	cm	inches
S	10.2 - 10.8	4" - 4.25"
M	10.8 - 11.4	4.25" - 4.5"
L	11.4 - 12	4.5" - 4.75"
XL	12 - 13.1	4.75" - 5"
XXL	13.1 - 14.1	5" - 6.25"

Positioning of Leatt® protection range

No single size of knee brace can accommodate all body dimensions because of variation in height and shape. In addition, the Leatt® Knee Brace X-Frame can only provide the maximum achievable level of protection when it is correctly positioned and adequately secured on the wearer. Where fitted the X-Frame should always be of a snug fit to the leg. In particular users should avoid choosing too big a protector because it can result in a dangerous riding condition. If the X-Frame does not fit correctly or if the knee protector is not positioned centrally on the knee cap and along the line of your shin such that correctly covers the knee area, an alternative size may be required.

Adjustment

Ensure that the X-Frame hinge is positioned in line with the centre of the knee on the side (in line with the "hinge" of the knee). Unfasten the adjustment straps. Pull the straps and adjust to a snug fit. Once the tension of the straps is adjusted, clips can be used to attach and release the straps during subsequent uses. Straps should be secured in the following order to ensure optimal fit:

When correctly-adjusted and worn the X-Frame should allow the full range of normal leg movement. With your knee brace X-Frame in place, check that you are able to safely operate your vehicle's controls. Please ask the retailer or contact Leatt® Corporation if you require further advice on sizing, adjustment or fitting.

Marking and labeling


The following diagrams provide important information about this product and the protection it provides.

CE	The "CE " mark denotes compliance with the health and safety requirements of Annex II of the Personal Protective Equipment Regulation (EU) 2016/425.
	The recommendation to read the manufacturer's instructions contained in this booklet
K TYPE A 1 EN 1621-1:2012	Denotes that the product is or incorporates a knee protector ("K - Type A") which meets level 1 of European Standard EN 1621-1:2012 The standard provides for two sizes of protector Type A and Type B representing minimum zones of protection and two levels of protection (impact attenuation), Level 1 and Level 2. Type B protector dimensions are intended to cover the needs of most motorcycle riders. However, for reasons of size, fit, ergonomics or comfort alternative (smaller) Type A protector dimensions may be used. Level of protection is determined by measuring the impact force transmitted in three defined strike areas of the minimum zone of protection, A B and C, after being subjected to a 50 Joule impact under ambient and wet conditions (optionally also hot and/or cold). Level 1 requires a mean impact force transmission of ≤ 35kN, all single strike forces to be ≤ 35kN in strike area A and ≤ 50kN in strike areas B and C. Level 2 requires a mean impact force transmission of ≤ 20kN,all single strike forces to be ≤ 20kN in strike area A and ≤ 30kN in strike areas B and C.

Cleaning

Wipe only with a soft damp cloth in gentle soap or mild detergent in cold water, rinse by wiping with a soft cloth damped in cold water and dry with a soft cloth. Comfort padding can be hand washed and air dried.

Common substances applied to this product (solvent, cleaners, hair tonics, etc.) can cause damage that may be invisible to the user and compromise the effectiveness/safety of your Leatt® product. If in any doubt, please contact the manufacturer before cleaning or laundering your protector. Always refer to the material content label on the Leatt® product. None of the materials used in the construction of the product are known to contain anything that is harmful or might cause an allergic reaction or be dangerous to health.

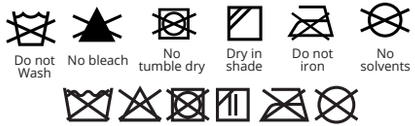Storage and transportation

Place the Leatt® product on a clean, dry surface and store in a dry, well-ventilated environment away from direct sunlight and extremes of temperature. Do not place heavy objects on top of Leatt® product as this can damage the fabrics and plastic materials. The protection provided by your Leatt® product should not be affected by non-extreme ambient temperature and weather conditions

Maintenance and obsolescence

For maximum performance, the Leatt® product must be inspected prior to each use. Stop using the protector if any components show signs of damage or are worn, cracked, and/or deformed, or if the inside padding has deteriorated. Leatt® Corporation recommends that the Leatt® product is replaced every three (3) years under normal use conditions. Should you experience a collision or if the Leatt® product or any component shows any signs of damage, cracking, deformity and/or the inside padding is deteriorating, return the product to the seller for inspection or destroy and replace it. Please dispose of it in an environmentally-friendly manner. Do not drill holes, paint, cut, heat, bend, apply decals or stickers, or make any other modifications to the knee brace. Doing so may compromise the effectiveness and safety of your Leatt® product.

Use with other Personal Protective Equipment

Leatt® protectors form part of a protective equipment ensemble available for use by participants in off road motorcycling and downhill biking activities. Additional PPE, which it is recommended should be worn during off road motorcycling and downhill mtb activities, includes clothing conforming to EN 13595-1:2002, where applicable to the riding environment (may not be suitable for off-road motorcycling or downhill mtb) or other suitable clothing providing full coverage to the body, arms and legs; impact protectors conforming to EN 1621-1:2012 (where not already a component of your Leatt® protector or other motorcycle clothing); gloves conforming to EN 13594:2002; boots conforming to EN 13634:2010 Always perform any motor sport activity with caution and use/wear appropriate safety equipment.

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Copyright © 2025, Leatt® Corporation
Leatt® Corporation HQ
Leatt® Corporation USA